

Calculate chlorophyll a concentration and Symbiodiniaceae counts

If the metadata sheet is filled out correctly, R functions, stored in a package called `coralchlo`, can be used to calculate the chlorophyll a concentration with the absorption measured with the photometer. The coral fragment area is calculated with the change in weight before and after adding the second layer of paraffin wax and used to normalize the chlorophyll a concentration and Symbiodiniaceae counts per area.

The package can be installed with the `devtools` package:

```
devtools::install_github("andieich/coralchlo")
```

The current version is 0.0.0.3.

After installation, you can load the package and download the metadata sheet. During this the import of the sheet, some basic test are done to ensure it was filled out correctly and that the values make sense.

Additionally, the surface area is estimated from the increase in weight with the additional layer of paraffin wax. Different conversion factors can be used. Either a value from the literature from Veal et al. (2010), or conversion factors based on custom-made calibration blocks (see [here](#) for more info). The method can be selected in the `read_metadata()` function with `method = "veal"` or `method = "criobe"`. The default selection is `"criobe"`.

```
library(coralchlo)

link_metadata_sheet <- "example_link"

dat_overview <- read_metadata(link_metadata_sheet)
```


Then, the chlorophyll a concentration is calculated based on the absorption measured with the photometer.

The photometer files (.csv) can be store in Google Drive or locally. If they are stored locally, `path_to_biotekfolder`, is the path to the folder containing these files. If they are stored in Google Drive, `path_to_biotekfolder` is the path within the Google Drive to the folder containing the photometer files. If you use Google Drive, set `is_googledrive = TRUE`. The files are downloaded to a temporary folder or to a folder specified in `download_directory`. All files in `download_directory` will be replaced.

The values are blank-corrected. To ensure that the blank values make sense, they are plotted for all wavelengths and photometer files. This can be skipped by setting `plot = FALSE` in the `normalise_chl_per_area()` function. The area data is used to normalize the chlorophyll concentration for the surface area of the coral fragment. For each sample, two measurements were taken, therefore two chlorophyll a and c2 concentrations will be exported.

sample_id	measurement_replicate	chl_a_per_cm2	chl_c2_per_cm2	chl_tot_per_cm2
2TSML27P	m1	2.514030	0.3057623	2.819792
2TSML27P	m2	2.514030	0.3057623	2.819792
2TSPL07P	m1	2.013116	0.2365427	2.249658
2TSPL07P	m2	2.017770	0.2350646	2.252834
2TSAL09P	m1	3.948941	0.4135201	4.362461
2TSAL09P	m2	3.948941	0.4135201	4.362461
2TSML26P	m1	4.847892	0.4351625	5.283054
2TSML26P	m2	4.847892	0.4351625	5.283054
2TSAL07P	m1	4.706739	0.4407329	5.147472
2TSAL07P	m2	4.706417	0.4543703	5.160787

```
dat_chl <- normalise_chl_per_area(dat_overview,
                                  path_to_biotekfolder = "example_path")
```

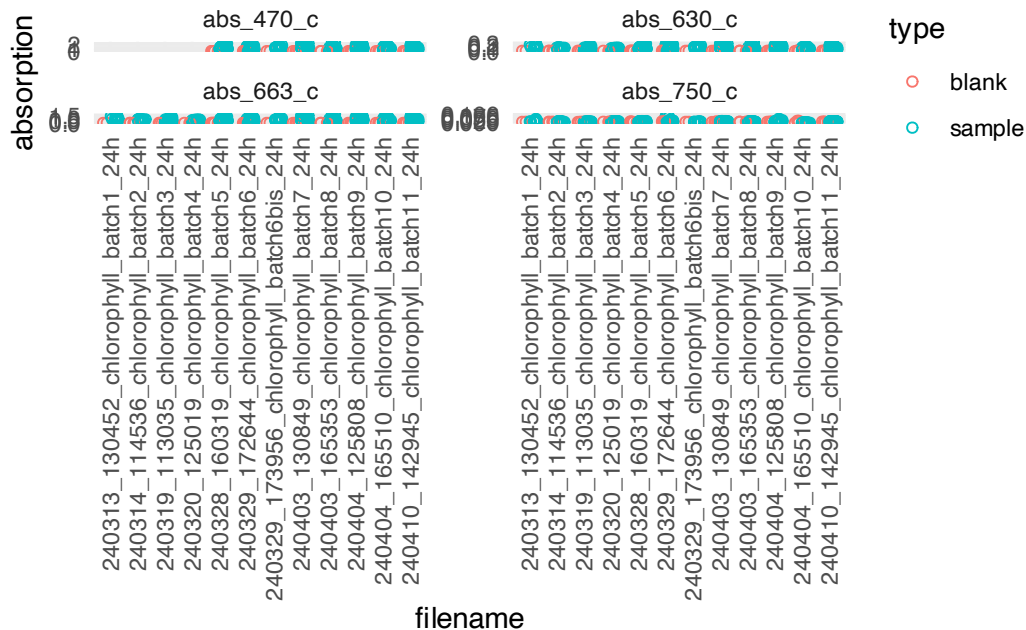

```
head(dat_chl, n = 10)
```

Similarly, the Symbiodiniaceae counts are normalized for the surface area of the coral fragment. Since six measurements are taken per sample, six values are exported

```
dat_counts <- normalise_counts_per_area(dat_overview)
head(dat_counts, n = 10)
```

sample_id	count_replicate	count_per_cm2
2TSML27P	c1	3623635
2TSML27P	c2	3355218
2TSML27P	c3	4026261
2TSML27P	c4	3221009
2TSML27P	c5	4160470
2TSML27P	c6	3489426
2TSPL07P	c1	3868471
2TSPL07P	c2	4479282
2TSPL07P	c3	3766669
2TSPL07P	c4	3155858

Veal, C. J., Holmes, G., Nunez, M., Hoegh-Guldberg, O., & Osborn, J. (2010). A comparative study of methods for surface area and three-dimensional shape measurement of coral skeletons. *Limnology and Oceanography: Methods*, 8(5), 241–253. <https://doi.org/10.4319/lom.2010.8.241>